

食品分析

3-2

12 种 C4~C18 饱和/不饱和脂肪酸甲酯的分析

组份: C4、C6、C7、C8、C10、C12、C14、C16、C18:0、C18:1、C18:2、C18:3 脂肪酸甲酯


色谱柱: MEGA-WAX 25m x 0.32mm 0.25um df

订货号: S-WAX-032-025-25



FAMES C4 - C18:3 (Fatty Acids Methyl Esters)

Column: **MEGA-WAX** - 0.32mm, 0.25µm, 25m
Catalog Code: S-WAX-032-025-25

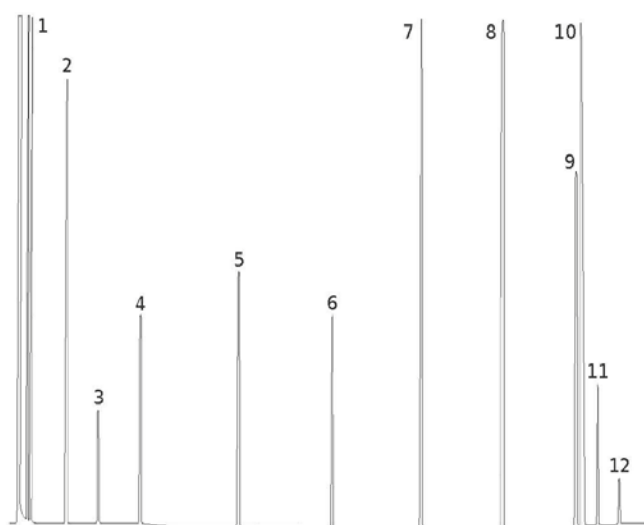

Conditions:

Injection: Split 250°C, 1.0µL, 1:50 Split Ratio.

Detector: FID 250°C.

Oven Temp.: 50°C, 4°C/min, 210°C.

Carrier Gas: Hydrogen, 60 kPa.

Peaks:

1. Butyric Acid ME (C4)
2. Caproic Acid ME (C6)
3. Heptanoic Acid ME (C7)
4. Caprylic Acid ME (C8)
5. Decanoic Acid ME (C10)
6. Lauric Acid ME (C12)
7. Myristic Acid ME (C14)
8. Palmitic Acid ME (C16)
9. Stearic Acid ME (C18:0)
10. Oleic Acid ME (C18:1)
11. Linoleic Acid ME (C18:2)
12. Linolenic Acid ME (C18:3)

since
1980

MEGA s.n.c. - Capillary Columns Laboratory - Via Plinio, 29 - 20025 Legnano (MI) - Italy - Tel./Fax: (+39) 0331.54.79.24 - info@mega.mi.it

www.mega.mi.it

