



农药分析

2-2

## 7 种有机磷农药的分析

组份：二嗪农、甲基对硫磷、甲基毒死蜱、马拉硫磷、乙硫磷

色谱柱：MEGA-1701 25m x 0.25mm 0.30um df

订货号：S-1701-025-030-25



### Organophosphorus Pesticides on MEGA-1 & MEGA-1701

Column: **MEGA-1701 - 0.25mm, 0.30 $\mu$ m, 25m**

Catalog Code: S-1701-025-030-25

#### Conditions:

Injection: Split 250°C, 1.0 $\mu$ L, 1:20 Split Ratio.


Detector: FPD 160°C.

Oven Temp.: 160°C (1min), 5°C/min, 250°C.

Carrier Gas: Hydrogen, 60kPa.

#### Peaks:

1. Diazinone
2. Chlorpyrifos methyl
3. Parathion methyl
4. Malathion
5. Methidathion
6. Ethion
7. Azinphos methyl

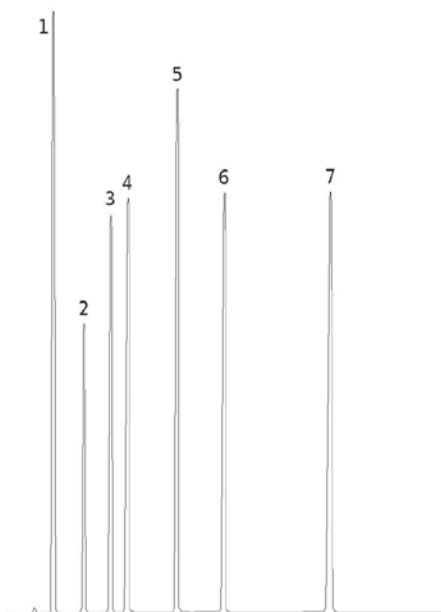

Acknowledgement: Prof. C. Bicchì, Dr. A. D'Amato et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.

since  
1980

MEGA s.n.c. - Capillary Columns Laboratory - Via Plinio, 29 - 20025 Legnano (MI) - Italy - Tel./Fax: (+39) 0331.54.79.24 - info@mega.mi.it

www.mega.mi.it

